

Freedom to Breathe

Aero

PISTON COMPRESSOR NEBULIZER



- Light Weight
- Innovative Design
- Compact - Easy To Carry
- Made From High Impact Strong & Durable Plastic

AEROSOL THERAPY UNI

Aero


PISTON COMPRESSOR NEBULIZER

Aero nebulizer is designed for continuous use built with heavy duty oil free piston compressor. No Maintenance is required as the oil free piston compressor does not need any lubrication.

Through the **Aero** nebulizer, a treatment of 2cc of medicine is completed in 7 minutes, nebulized in very small particles measuring between 0.5 and 5.5 microns, reaching the deepest areas of the respiratory tract.

Dust and impurities are kept away by the filter which ensures a safe medical use & is easily replaceable.

Each **Aero** nebulizer is thoroughly tested and is covered with 1 year warranty.

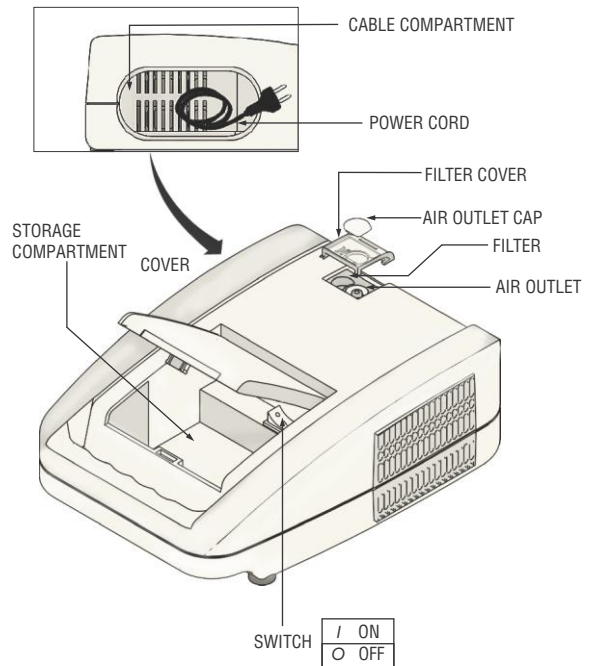

WARRANTY:

Each **Aero** nebulizer is carefully inspected and is covered with 12 months warranty.


THE **Aero** NEBULIZER IS SUPPLIED WITH:

- Nebulizer Cup
- Mouthpiece
- Mask
- Air filter
- Medical grade Nebulizer tubing

SPECIFICATION:

Electricity Supply	: 230 V ~ 50 Hz
Pressure	: 1.5 bar
Airflow	: 8 lt./min.
Size (mm)	: 207 x 273 x 97
Weight	: 2.2 Kg. approx

* The company reserves the right to make technical and functional changes and other accessories without prior notice.

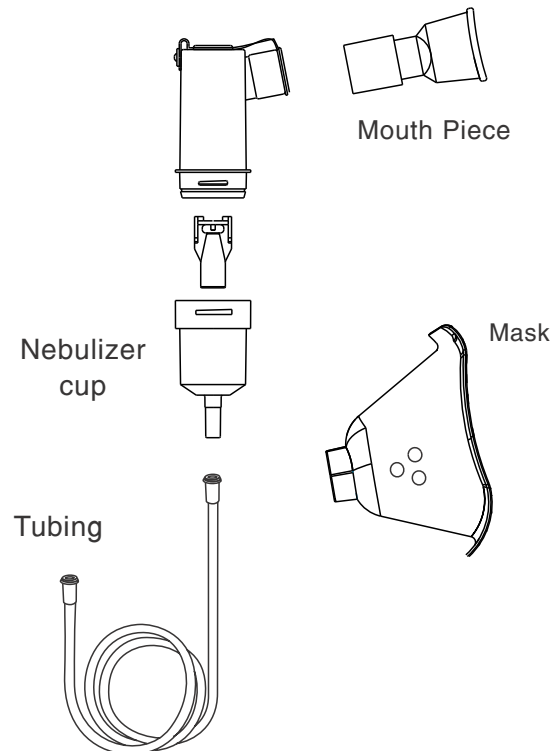

Manufactured & Marketed By :

HEMANT SURGICAL INDUSTRIES LTD.

647, Girgaum Road, Avijah Building, Shop No. 6 & 7, Near Metro Cinema, Mumbai-400 002. • Tel. : 2203 0935, 2203 0936 • Fax : 91-22-2207 7585
Website : www.hemantsurgical.com • E-mail : info@hemantsurgical.com